



New! Forced Convection Lab Ovens

Programmable Forced Convection Lab Ovens

Yamato DKV400/410/600/610/810/910

MADE IN



For Information and to Order
Contact Us At:
323-770-0634 800-574-2748
sales@LRE.com www.LRE.com

LR
TECHNOLOGIES

Operating temp.
range

RT+10°C~270°C
DKV400/410/600/610

RT+10°C~220°C
DKV910

Temp. fluctuation

±.5°C - ±.6°C

Internal
capacity

3.1 CuFt
(DKV400/410)

5.3 CuFt
(DKV600/610)

10.6 CuFt
(DKV810)

19.2 CuFt
(DKV910)

Operation and functions

- Superior temperature accuracy
- Highly visible 3.5 inch color LCD controller (VSD-1 type)
- Door hinges with click function at 30°click and 120°click
- Fixed temperature operation, Quick Auto-Stop Operation / Auto-Stop Operation / Auto-Start Operation / Programmed Auto-Start Operation (Time and Clock Settings), Programmed Operation (Time Settings)
- Increased safety and self-diagnostic function
- Easy to use and maintain
- Built in exhaust ports

Q: Is it possible to install an observation window on the door?

A: Custom orders can be accommodated - however, insulation and specified performance may vary depending on operating temperature. Standard window size is 9.8" X 11". For alternative sizes, please let us know.



3.1 CuFt

DKV400 W/Optional Stand



5.3 CuFt

DKV600 W/Optional Stand

Specifications

Model	DKV400	DKV410	DKV600	DKV610	DKV810	DKV910
Circulation method	Forced convection					
Operating temp. range	Room temp. +10°C to 270°C					RT +10°C to 220°C
Temp. fluctuation ^{*2}	±.6°C (at 270°C)		±.5°C (at 270°C)		±.6°C (at 270°C)	±.6°C (at 220°C)
Temp. gradient ^{*2}	9°C (at 270°C)				11°C (at 270°C)	11°C (at 220°C)
Temp. rise rate ^{*2}	3.3°C / min		2.1°C / min		3.4°C / min	2.5°C / min
Temp. rise time ^{*2}	~75 min.		~115 min.		~75 min.	~80 min.
Exterior material	Chromium-free electrogalvanized steel sheet, baked-on finish					
Chamber	Stainless steel					
SUS tube heater	1.5 KW		2.5 KW		3.0 KW	
Cable port	Inner diameter 1.29" (1 piece on right side)					
Exhaust port	2 Pcs (1.3" I.D) Top				2 Pcs (1.3" I.D) Back	
Safety devices	Self-diagnostic functions (automatic overtemperature protection, temperature sensor error, SSR short circuit, heater disconnection, fan error, clock element error, main relay error, memory error, internal communication error) / independent overheat prevention device / door switch (door alarm) / overcurrent leakage breaker / fall prevention fitting (Model 810/910)					
Controller type	VSD-1 Controller					
Temp. control system	PID control by microcomputer					
Temp. display	3.5-inch color LCD panel digital display					
Timer / clock timer	1 minute to 99 hours 59 minutes / Date: hour: minute (up to 31 days later)					
Operation modes	Fixed temperature / Timer: Hourly quick auto stop operation / Timer: Hourly auto stop operation / Timer: Hourly auto start operation / Program operation / Timer: Hourly program auto start operation					
Additional functions	Calibration offset / Power failure return mode selection function / Timer Wait function / Security management function / Temperature upper/lower limit alarm function / Fan OFF delay function / Power failure alarm / External communication terminal (RS-485) / Temperature output terminal (4-20mA) / Time-up output terminal / Error output terminal / Recorder terminal / Program preview function / Operation trend function					
Sensor	Temperature controller / independent overheat prevention device: K thermocouple					
Safety devices ^{*3}	Self diagnosis functions (temp. sensor abnormal, Heater disconnection, SSR- short, Automatic overheating prevention), Key lock function, Overheat prevention, Electric leakage breaker with over current protection					
Internal dimensions (W×D×H)	17.7" x 17.7" x 17.7"		23.6" x 19.6" x 19.6"		23.6" x 19.6" x 39.4"	42.9" x 19.6" x 39.4"
External dimensions (W×D×H) ^{*4}	22.8" x 25" x 32.7"		28.7" x 27" x 34.6"		28.7" x 27" x 63.3"	48" x 27" x 63.5"
Internal capacity	3.1 CuFt		5 CuFt		10.6 CuFt	19.2 CuFt
Chamber shelf load	33 Lbs /Shelf					
Chamber load capacity	66 Lbs.		66 Lbs.		132 Lbs	264 Lbs.
Number of shelves / pitch	11 steps / 1.1"		13 steps / 1.1"		29 steps / 1.1"	29 / 1.1" x2 columns
Power supply 50/60Hz Rated current	115V, 13.5A with plug	220V, 7.5A no plug, round terminal	115V, 13.5A with plug	220V, 7.5A no plug, round terminal	Single phase 220V, 12A with plug	Single phase 220V, 14.5A no plug, round terminal
Weight	132 Lbs		165 Lbs		264 Lbs	429 Lbs
Included shelf plate / brackets	1 pc. / 2 pcs.				3 pc. / 6 pcs.	6 pc. / 12 pcs.
Fall prevention fittings	-		-		4 pcs (fixing hole diameter Φ14 M12 recommended)	

^{*1} Performance (excluding temperature rise time) is measured at rated voltage, 50Hz or 60Hz frequency, room temperature 23° C ±5° C, humidity 65%RH ±20%, and no load. Note that the maximum temperature in the temperature control range may not be reached if the power supply voltage is below the rating or if the outside temperature is low. ^{*2} Temperature rise time is reference data for rated voltage, 50 Hz or 60 Hz, room temperature 23° C, humidity 65%, and no load. Note that the temperature rise time may be longer when the power supply voltage is below the rating or when the outside temperature is low. ^{*3} If the product is used beyond the operating voltage range (±10% of rating), there is a possibility that the safety function will malfunction. Pay attention to the power supply voltage when using the product. ^{*4} Protrusions are not included.



10.6 CuFt

DKV810




19.2 CuFt

DKV910

Control

Panel

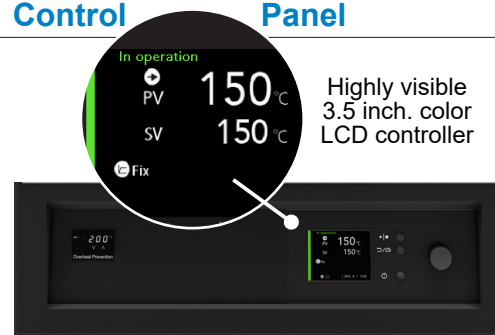

Interior



DKV400/410



DKV600/610



DKV810





DKV910

Standard DKV910 model includes a center column, as shown in the photo. However, we can accommodate special orders for versions without the center column.

*Please note that the maximum temperature will be limited.

Door Hinges With Click Function

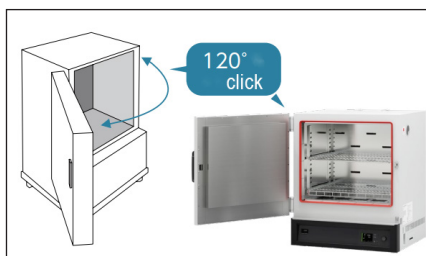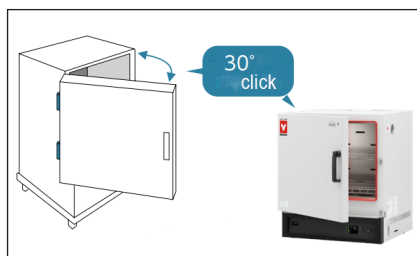

Door hinges have a click function that can hold the door at two angles: 30° and 120° (operating range is 180°). Samples can be inserted and removed safely, without concern for burns or accidental contact.

Optional Items

Description	Factory Code	Part Number	Compatible Models
Stand (no casters)	OA194	281596	DKV400/410/600/610
	OA117	281355	DKV400/410
	OA118	281356	DKV600/610
Stacking kit	OA196	281598	DKV400/410
	OA197	281599	DKV600/610
	-	212246	DKV400/410
Shelf and bracket set Stainless steel perforated metal shelf. Same shelf as standard accessories.	-	212266	DKV600/610/810
	-	212249	DKV910
	-	296902	DKV400/410/600/610
Seismic isolation mat Secures devices in place by fastening them to the base, preventing tipping			
Sheath sensor (K thermocouple)	ODT48	212946	DKV400/410/600/610/810/910
Silicone stopper Silicone rubber plug for cable hole, φ2mm hole for sensor	ODT52	212947	DKV400/410/600/610/810/910
*Rubber feet Rubber feet attachable to the main unit legs	OA198	281634	DKV400/410/600/610
*Cable port Contact Customer Service for installation position and quantity. Not available for installation after delivery.			
*Cable port 25mm inner Ø	ODM36	281454	DKV400/410/600/610/810/910
*Cable port 50mm inner Ø	ODM38	281455	DKV400/410/600/610/810/910



- Never use in flammable or explosive gas atmosphere.
- Never use explosive or flammable material.

- Caution: High temperature components.